


# Capture with Precision

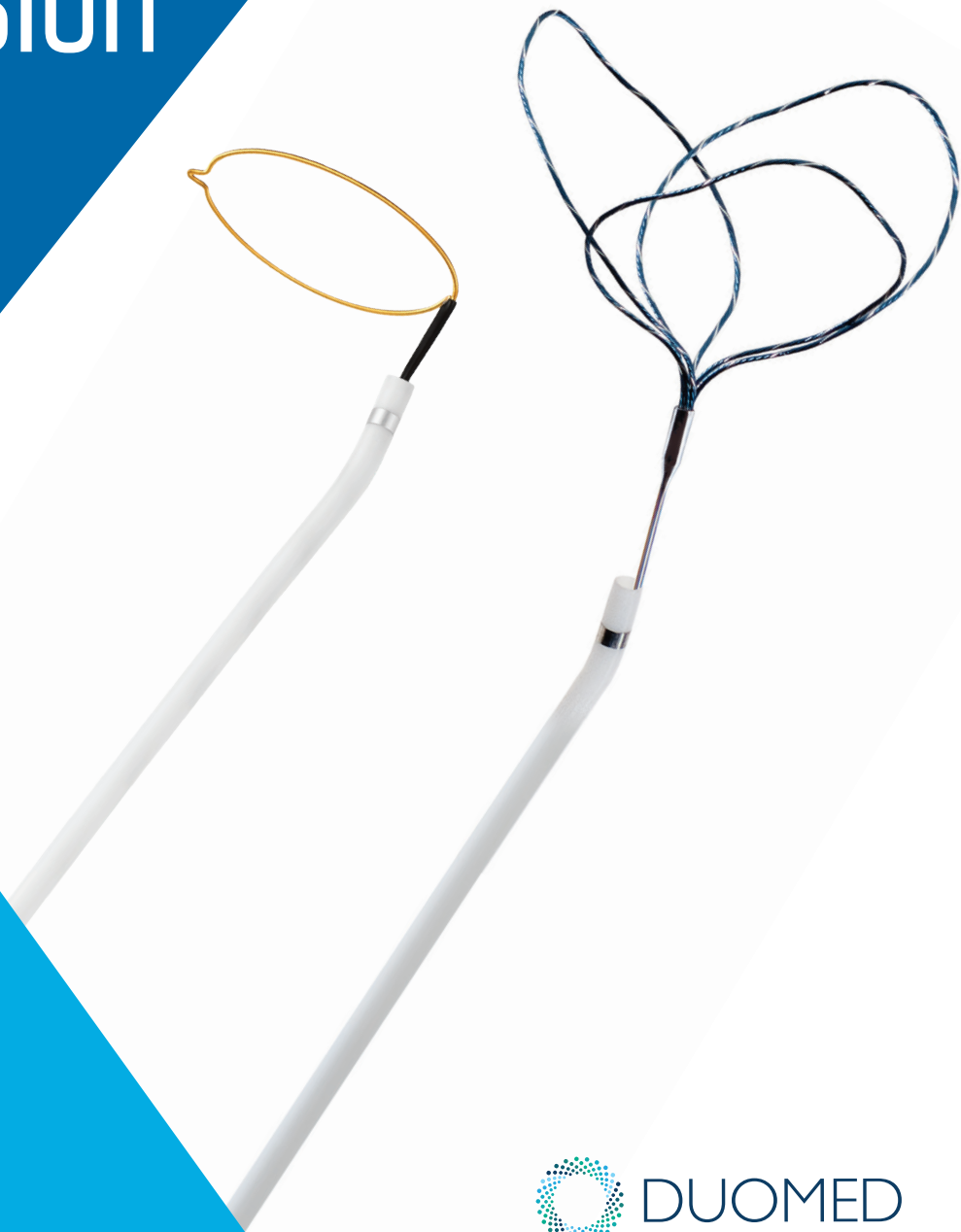

**HALO**<sup>™</sup>  
SINGLE-LOOP SNARE KIT

**Atrieve**<sup>®</sup>  
VASCULAR SNARE KIT

 DUOMED

  
**ARGON**  
MEDICAL DEVICES



The Atrieve Vascular Snare Kits are designed with 3 independent loops to provide a distinct advantage in the manipulation and retrieval of foreign bodies in the cardiovascular system or hollow viscus.

**Independent Nitinol Loops**

designed with short profile, tipped loops to improve capture through complete vessel coverage

**Heightened Visualization**

provided by the platinum strands incorporated into each loop

**Radiopaque Marker Band**

identifies the snare delivery catheter's location under imaging


**Improved Directional Control**

enabled by the 15° curved tip on 6F and 7F snare delivery catheters

**Torque Device**


Easy to grip design that provides 1:1 response when manipulating the snare during the procedure

\*Featured on Atrieve and Halo products

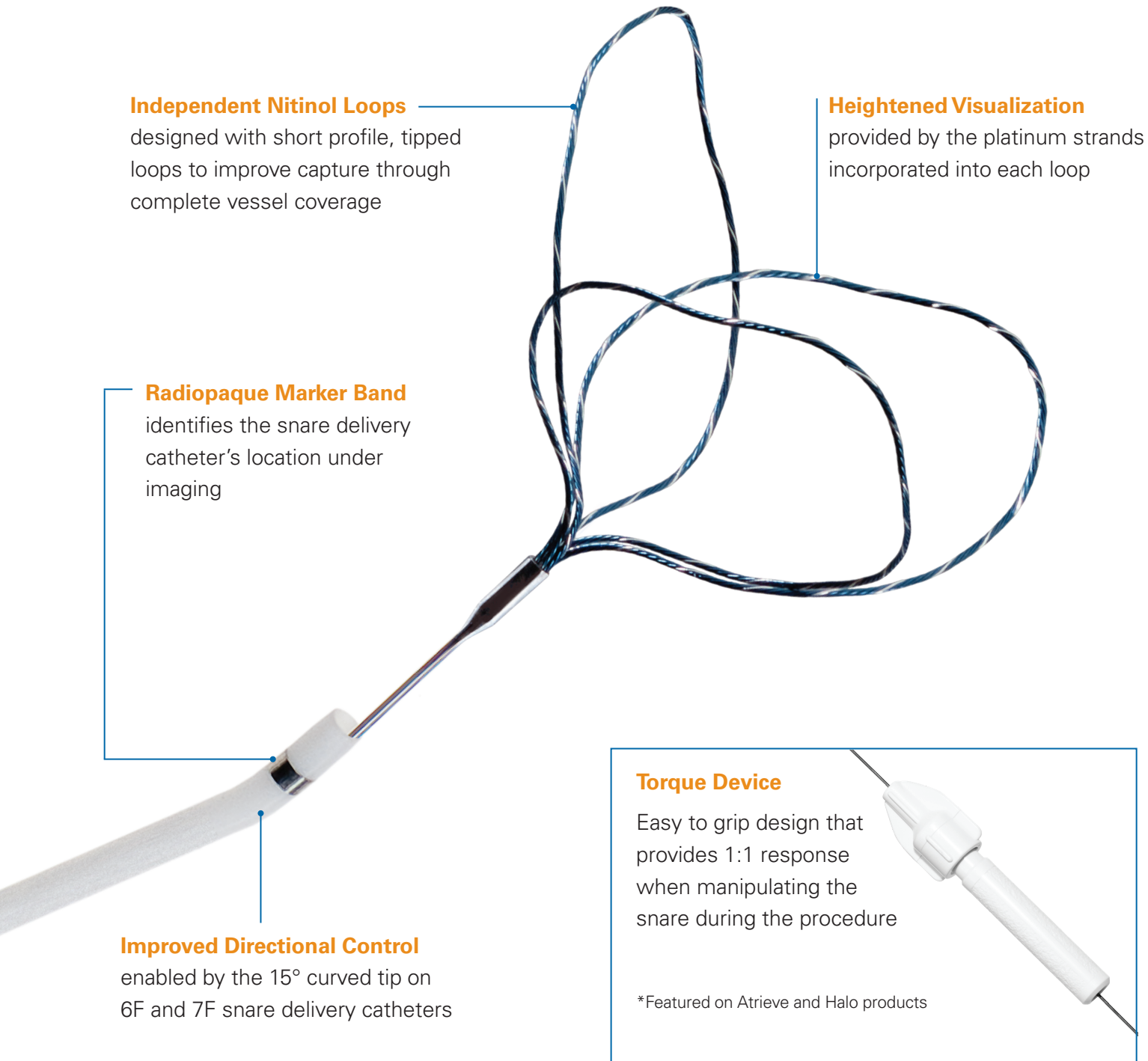




The Halo Single-Loop Snare Kits are designed to effectively retrieve and manipulate objects in the cardiovascular system or hollow viscus. Their 90° angled loop and nitinol design, allow for optimal wall apposition and ideal capturability.

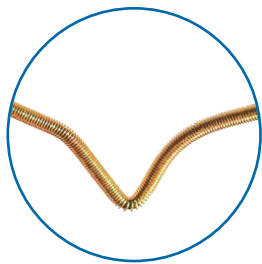

**Nitinol & Gold-Plated Tungsten Loop**

retains snare's circular shape, while allowing for superior visibility under fluoroscopy

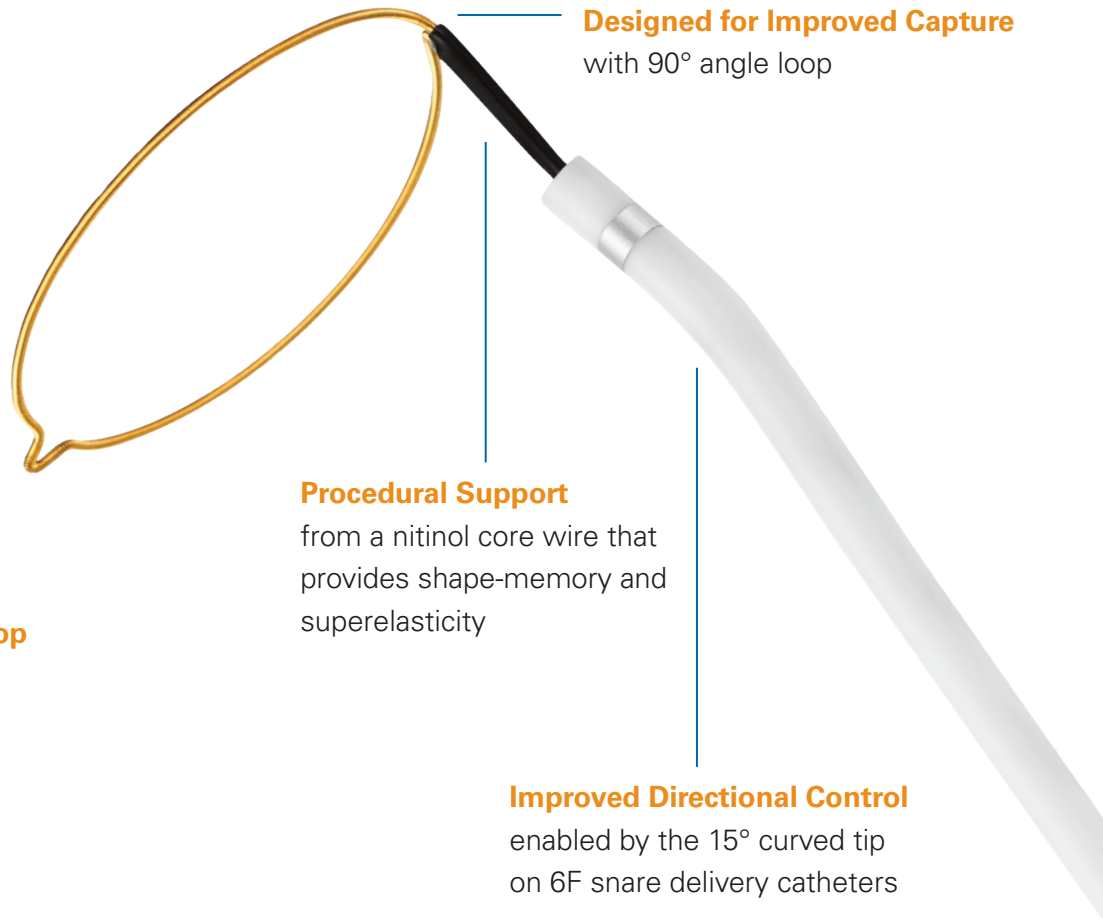

**Designed for Improved Capture**  
with 90° angle loop

**Procedural Support**

from a nitinol core wire that provides shape-memory and superelasticity

**Improved Directional Control**

enabled by the 15° curved tip on 6F snare delivery catheters

**Efficient Procedural Assembly**

An insertion tool with a unique reusable, peel away design that allows for quick and accurate snare loading



## Ordering Information

Atrieve® Micro Vascular Snare Kit (Box of 1)				
CATALOG #	LOOP DIAMETER	SNARE LENGTH	CATHETER SIZE	CATHETER LENGTH
381003004	2-4 mm	175 cm	3.2F	150 cm
381003008	4-8 mm	175 cm	3.2F	150 cm

Atrieve® Standard Vascular Snare Kit (Box of 1)				
CATALOG #	LOOP DIAMETER	SNARE LENGTH	CATHETER SIZE	CATHETER LENGTH
382006010	6-10 mm	120 cm	6F	100 cm
382006015	9-15 mm	120 cm	6F	100 cm
382006020	12-20 mm	120 cm	6F	100 cm
382007030	18-30 mm	120 cm	7F	100 cm
382007045	27-45 mm	120 cm	7F	100 cm

Halo™ Single-Loop Microsnare Kit (Box of 1)				
CATALOG #	LOOP DIAMETER	SNARE LENGTH	CATHETER SIZE	CATHETER LENGTH
383302175	2 mm	175 cm	3F	150 cm
383302200	2 mm	200 cm	3F	175 cm
383304175	4 mm	175 cm	3F	150 cm
383304200	4 mm	200 cm	3F	175 cm
383307175	7 mm	175 cm	3F	150 cm
383307200	7 mm	200 cm	3F	175 cm

Halo™ Single-Loop Snare Kit (Box of 1)				
CATALOG #	LOOP DIAMETER	SNARE LENGTH	CATHETER SIZE	CATHETER LENGTH
383405120	5 mm	120 cm	4F	102 cm
383410065	10 mm	65 cm	4F	48 cm
383410120	10 mm	120 cm	4F	102 cm
383615120	15 mm	120 cm	6F	102 cm
383620120	20 mm	120 cm	6F	102 cm
383625065	25 mm	65 cm	6F	48 cm
383625120	25 mm	120 cm	6F	102 cm
383630120	30 mm	120 cm	6F	102 cm
383635120	35 mm	120 cm	6F	102 cm



**DUOMED**

**Duomed Swiss AG**  
 Grenzstrasse 5a, 6214 Schenkon  
 Tel. +41 (0)41 510 07 00  
[www.duomed.com](http://www.duomed.com)

1. K021606, K200268

© 2021 Argon Medical Devices, Inc. All rights reserved. Argon Medical Devices, the stylized logo, and the product and service marks mentioned herein are trademarks or registered trademarks of Argon Medical Devices, Inc. in the United States and other countries. All other trademarks are the property of their respective owners.

2600 Dallas Parkway, Suite 440, Frisco, TX 75034  
 P: 800.927.4669 F: 903.677.9396 [argonmedical.com](http://argonmedical.com)



96-6300-12/102021F