



DUOMED

Duomed Swiss AG
Grenzstrasse 5a, 6214 Schenkon
Tel. +41 (0)41 510 07 00
www.duomed.com

Thrombus Management

- **Option® ELITE Retrievable Vena Cava Filter**
- **CLEANER XT® and CLEANER 15® Rotational Thrombectomy Systems**
- **Atrieve® Vascular Snare Kit**




Option® ELITE Retrievable Vena Cava Filter

Versatility & Proven Performance

The unique design of the Option ELITE IVC filter offers consistent and proven clinical performance. The filter can be deployed using a standard or over-the-wire technique from antecubital, jugular, femoral, or popliteal access points. Option ELITE can also be retrieved when the patient no longer requires PE protection.

Box of 1

Catalogue #	Description	Cava Diameter	Sheath Length	Sheath Size (ID)
352506070E	Femoral or Jugular Access	up to 30mm	70cm	5F
352506100E	Antecubital or Popliteal Access	up to 30mm	100cm	5F

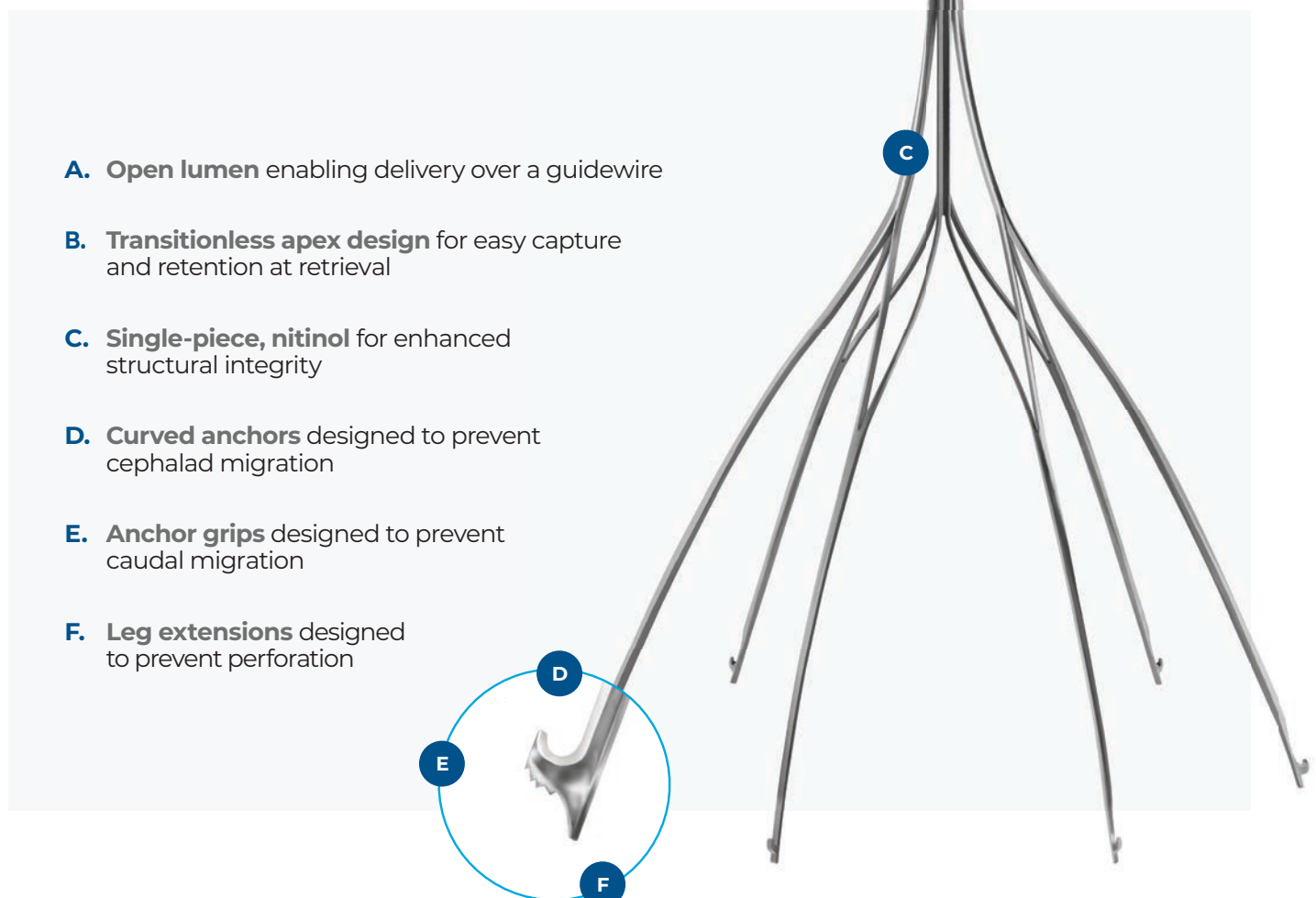


System Components

Curved Pusher

Deployment mark for alignment during delivery with a standard technique

Curved pusher aids centering and resists tilt during standard delivery

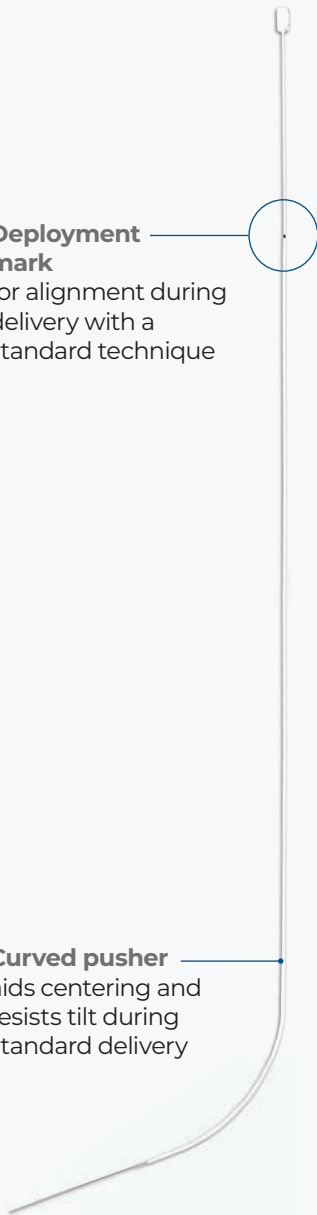

Delivery Sheath

Hemostatic valve to prevent bleed back

Integrated side port for hand injection

Universal Colour-coded filter cartridge for placement via four different access sites

Valve cap maintains hemostasis during post-placement cavagrams

Low-profile design with a 5 FR ID



Dilator

Deployment mark for alignment during delivery with an over-the-wire technique

Radiopaque markers spaced 32mm apart



CLEANER XT® and CLEANER 15® Rotational Thrombectomy Systems

The only wall-contacting rotational thrombectomy system indicated for mechanical declotting and infusion of physician-approved fluids for use in the peripheral vasculature, native vessel dialysis fistulae, and synthetic dialysis access grafts. CLEANER XT and CLEANER 15 Rotational Thrombectomy Systems are completely disposable and arrive in a single-piece system, available for use immediately.

CLEANER XT

For greater steerability and effective maceration in smaller lumen vessels




Box of 1

Catalogue #	Catheter OD	Length	Sinusoidal Wire Amplitude	Wire Diameter
700009XT	6F	65cm	9mm	.035"
700109XT	6F	135cm	9mm	.035"


CLEANER 15

For powerful and effective maceration in larger lumen vessels

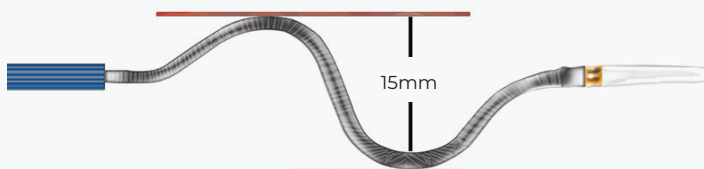

Box of 1

Catalogue #	Catheter OD	Length	Sinusoidal Wire Amplitude	Wire Diameter
700015	7F	65cm	15mm	.044"
700115	7F	135cm	15mm	.044"

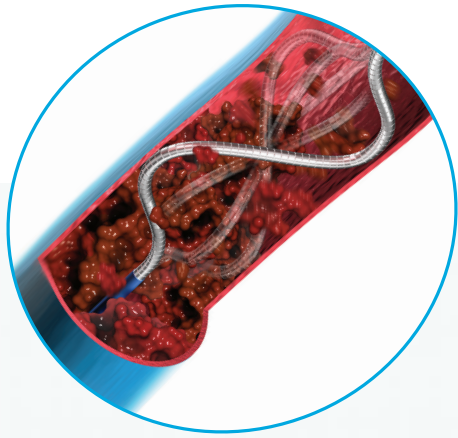


Sliding lever
collapses the
sinusoidal
wire for easy
advancement



**Maceration
button**
activates the
sinusoidal wire


3-way side infusion
allows infusion of
thrombolytics or
contrast media

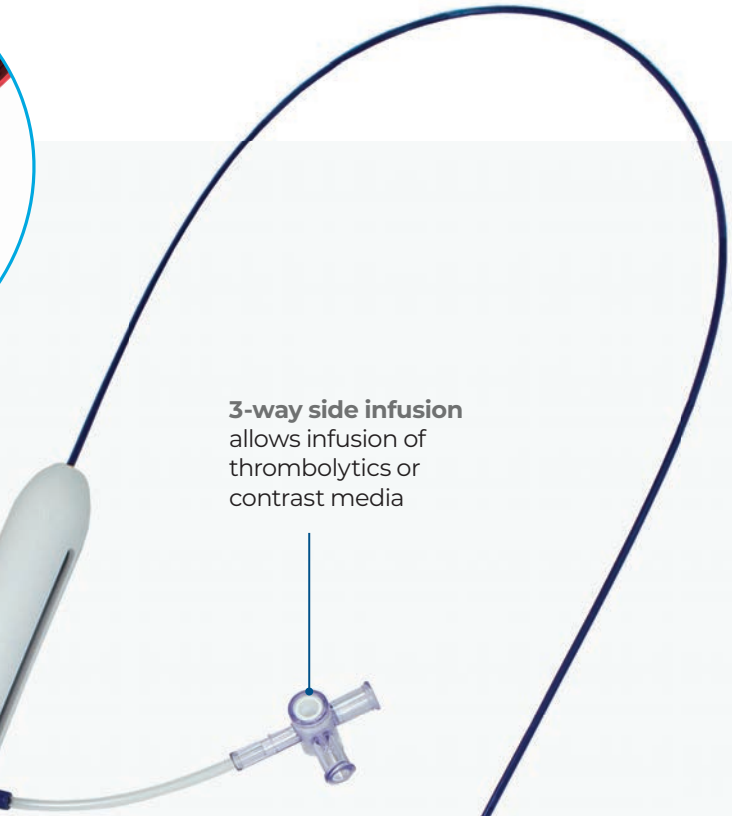

Distal side hole
delivers fluids
infused through
the side port


**Self-sizing
sinusoidal wire**
macerates thrombus
and restores vessel
patency (available
in 9mm and 15mm
amplitudes)

Atraumatic radiopaque tip
enables clear visibility
under fluoroscopy

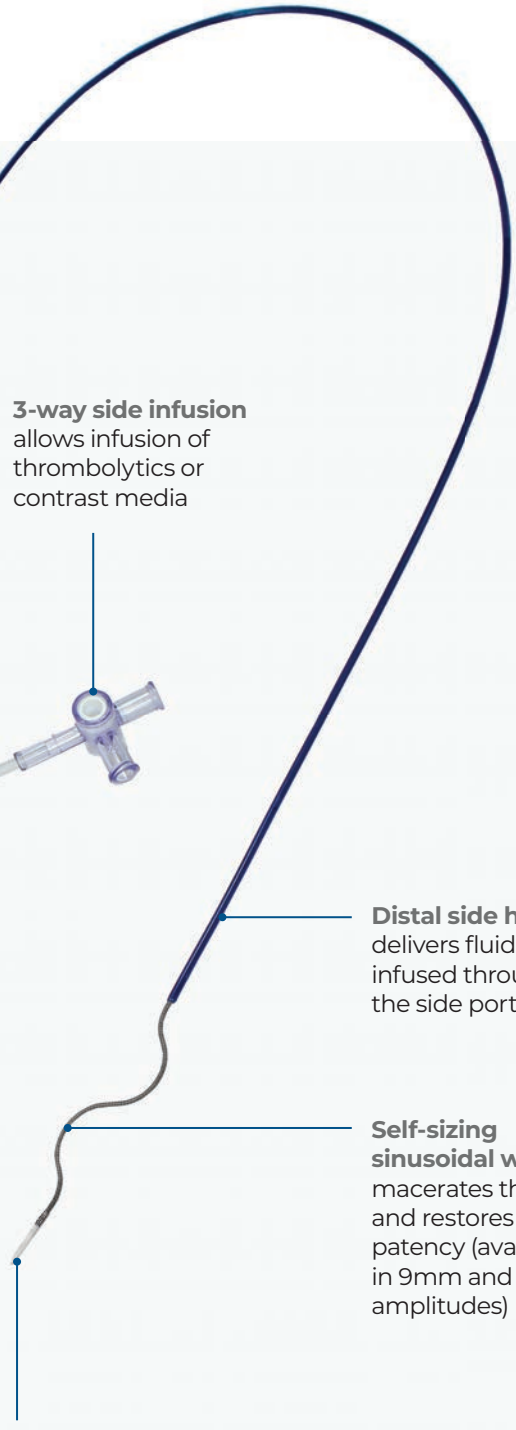

Atrieve® Vascular Snare Kit

The Atrieve Vascular Snare Kits are designed with three independent loops to provide a distinct advantage in the manipulation and retrieval of foreign bodies in the cardiovascular system or hollow viscus.

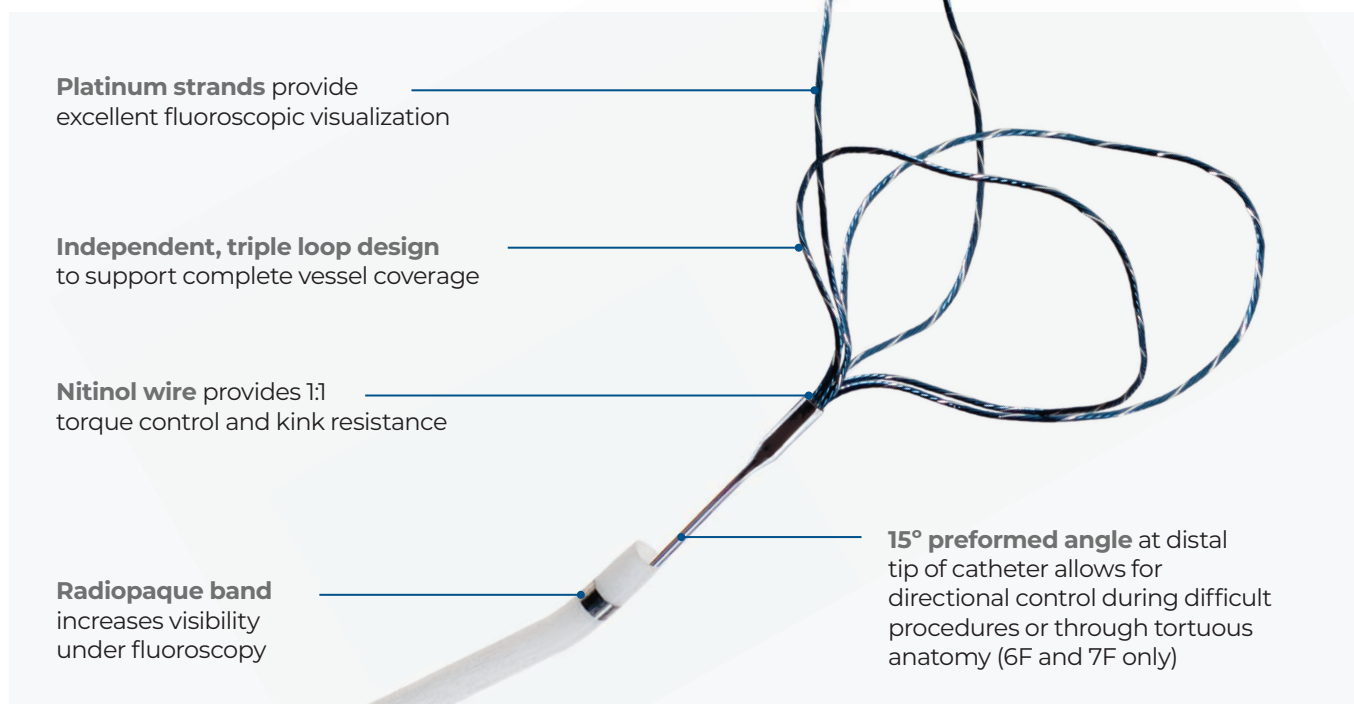

Box of 1

Catalogue #	Loop Diameter	Snare Length	Catheter Size	Catheter Length
381003004	2-4mm	175cm	3.2F	150cm
381003008	4-8mm	175cm	3.2F	150cm
382006010	6-10mm	120cm	6F	100cm
382006015	9-15mm	120cm	6F	100cm
382006020	12-20mm	120cm	6F	100cm
382007030	18-30mm	120cm	7F	100cm
382007045	27-45mm	120cm	7F	100cm



DUOMED

Duomed Swiss AG

Grenzstrasse 5a, 6214 Schenkon

Tel. +41 (0)41 510 07 00

www.duomed.com